COMPLETE SYSTEM User Manual

Deutsch
English
Français
Español
Italiano
Magyar

[Art. No.: 001211]

BEMER[®]
3000

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Operating instructions

Congratulations on acquiring the BEMER complete system. This high quality product was produced by us with the utmost care in accordance with our high quality standards and the European production standards that are applicable for medical products. In combination with user friendly functionality, our high quality standards have, for many years, given our customers above average satisfaction.

BEMER is good for you – stay healthy with BEMER.

Manufacturer:

INNOMED International AG

Schliessa 12 // FL-9495 Triesen

Made in Principality
of Liechtenstein.

CE 1275

In order to achieve the best possible therapy results, please observe the application options recommended for effective therapy (e.g. basic plan) that are described in the operating instructions.

User instructions // Foreword	50
Explanation of symbols	51
Use in accordance with intended purpose	52
Safety instructions	53
Contents	56
Legend	57
Start up procedure	60
Operation	62
User instructions	68
Functional test	79
Cleaning, servicing, disposal	81
Troubleshooting	84
Warranty	86
Technical data	88
Accessories	92

1. User instructions

These operating instructions contain important information. Keep these instructions with your unit so that you can refer to them at any time to check on safety instructions and correct handling of the unit.

Please read these operating instructions carefully before starting to use the unit. This will enable you to make use of all the advantages offered by the system and avoid any damage or injury to yourself and others.

The BEMER complete system is intended only for the purpose it has been designed as described in the operating instructions. Where these instructions have been translated into a foreign language the German original version shall be legally binding.

2. Foreword

BEMER electromagnetic field therapy is based on supporting fundamental physiological processes and hence on improving natural self-regulating mechanisms, as well as scientifically proven effects: effects on microcirculation, the immune system, protein biosynthesis and the vegetative nervous system.

Based on this alone, the BEMER therapy offers a very wide application spectrum with very good effects in terms of therapy, regeneration, rehabilitation and prevention.

The **BEMER *plus* signal** is based on the uniquely broad frequency spectrum of the BEMER impulse which has been used very successfully in therapy for many years. In addition to the already known physiological effects of the BEMER impulse, we have succeeded in activating the vasomotion of blood vessels significantly better than before through a certain modulation of the flux density (intensity) within the pulse sequence.

The term vasomotion refers to rhythmical changes in the diameter of arterioles and small arteries which play an important part in the distribution of blood in the tissue.

3. Explanation of symbols

Symbols in the operating instructions

Warning / Caution

Failure to observe this might cause risk of injury or damage to property.

Disposal

Note

Symbol on the unit

CE conformity marking

Agency reference 1275

Take note of user instructions

Protection rating

Unit with protection rating II

Type of application part

Application part of type BF

4. Use in accordance with intended purpose

The BEMER 3000 and BEMER 3000 *plus* are capable of stimulating the vasomotion (rhythmical contraction movements) of small and very small blood vessels of living organisms through weak pulsating electromagnetic fields generated by a special signal. This achieves a more favorable distribution of blood particles in the microcirculation of treated organisms, which in turn leads to an improved function of natural regulatory mechanisms, strengthening of the immune system, improved production of the body's own radical scavengers and an influence on protein synthesis. The BEMER 3000 systems can therefore be used in a preventative context as well as a com-

plementary therapy for the treatment of illnesses, in which the improvement of the microcirculation and any resulting effects on living organisms is therapeutically significant.

The BEMER complete system may not be used for any other purpose than that described in this section. Any use beyond this is considered not to be in accordance with the unit's intended purpose. The manufacturer accepts no liability for any injury or damage resulting there from. The risk rests solely with the user.

5. Safety instructions

Note

The BEMER complete system has been manufactured using state of the art technology and is safe to operate. Nevertheless, some risks remain. This is particularly the case if the unit is operated by insufficiently trained staff or if it is used improperly and not in accordance with its intended purpose.

- **Any conditions or illnesses must be assessed by a physician before using the BEMER complete system, should this not yet have been done.**
- **Any medication prescribed by a physician should not be changed unless confirmed by the treating physician.**
- **Any unexplained physical reaction occurring during a BEMER treatment should be discussed with the treating physician in order to be able to respond appropriately to an illness that may be developing.**
- **People with organ transplants, allergenic stem cell and bone marrow transplants are not permitted to use the BEMER complete system.**
- **People suffering from infectious diseases, severe cardiac arrhythmia, large aneurisms, severe psychoses and non-compensated seizure disorders (e.g. epilepsy) must obtain the approval of the treating physician before beginning treatment with the BEMER complete system.**
- **People with electronic implants should not use the BEMER complete system until it can be confirmed that the limits specified by the manufacturer of the implant are not exceeded.**

It is not possible to generally preclude any interactions with, and / or adverse effects on, electronic implants such as cardiac pacemakers, ICD's, brain pacemakers, insulin pumps etc.

- **Before connecting the BEMER complete system to the power supply, please ensure that the power supply voltage is identical to that stated on the power supply unit.**
- **The power supply cable and plug must be in good condition. Any damaged parts must be replaced immediately!**
- **Remove the plug from the power supply before carrying out any cleaning!**
- **Do not use the BEMER complete system for a purpose for which it is not intended.**
- **Always keep the BEMER complete system in proper order! Neglected units could lead to accidents.**
- **Regularly check the BEMER complete system for any damage. Any damaged parts must be replaced immediately!**
- **Take note of user instructions.**
- **The BEMER complete system may only be used on firm, hard surfaces (e.g. floor, couch, bed, table).**
- **The BEMER complete system is not a toy. Keep away from children.**

6. Contents

- 1 Towel
- 2 Operating instructions
- 3 BEMER Control (device to check the signal)
- 4 Intensive applicator
- 5 Power supply unit
- 6 Coil mat
- 7 Control unit
- 8 Pouch for intensive applicator
- 9 Carrying bag


7. Legend

7.1 Top of control unit

- 1 ON button
- 2 STOP button
- 3 Indicator light – polarity
- 4 Buttons for intensity levels 1 to 10
- 5 Buttons for programs P1 to P4
- 6 Connection socket for intensive applicator
- 7 Indicator light for intensive applicator
- 8 Time display for remaining therapy time
in 2-minute intervals from 20 to 0

7.2 Back of control unit

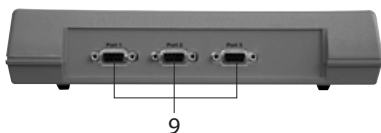- 9 Multifunction socket for power connection, coil mats or optional accessories, e.g. the Accu-Pak


7.3 Power supply unit

- 10 Power supply unit

- 11 Plug of power supply unit


7.4 Coil mat

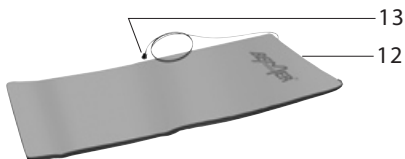- 12 Coil mat

- 13 Plug of coil mat

7.5 Intensive applicator

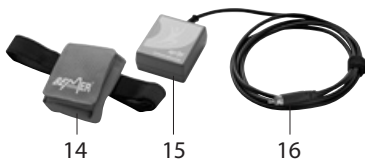- 14 Pouch for intensive applicator
- 15 Intensive applicator
- 16 Plug of intensive applicator

7.6 BEMER Control

- 17 Sensor range
- 18 Slide switch 'ON/OFF'
- 19 Control display

8. Start up procedure

You should regularly check your BEMER complete system, at least each time before starting the unit:

- The cable
- The control unit
- The power supply unit
- The coil mat
- The intensive applicator

The BEMER complete system may only be used or stored in rooms where it is protected from dust, humidity, dripping or splashing water.

The unit should be operated under the following ambient conditions:

- **Temperature: 10 °C to 40 °C / 50 °F to 104 °F**
- **Relative humidity: 30 to 75 % non-condensing.**


8.1 Control unit – coil mat

1. Connect the power supply unit plug to any socket on the control unit and attach firmly.

2. Connect the power supply unit to the electrical outlet.

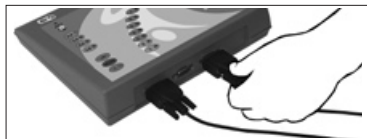3. Connect the coil mat plug to any socket on the control unit and attach firmly.

8.2 Control unit – intensive applicator

1. Connect the power supply unit plug to the control unit and attach firmly.

2. Connect the power supply unit to the electrical outlet.

3. Connect the intensive applicator plug to the control unit.

9. Operation

Note


The **BEMER *plus* signal** is the proven broad spectrum BEMER impulse with the same flux densities as before. Every 20 seconds, five pulses are emitted with a maximum flux density of approx. 35 microtesla (corresponds to level 10). This pulse can only be accessed through ports 1 to 3, using the programs P1 to P4. The coil mat and the intensive applicator cannot be used simultaneously. Although it is possible to connect them both at the same time, the coil mat will not function when the intensive applicator is plugged in. Applicators connected to the socket (6) have priority over the coil mat which is connected to the back.

9.1 Control unit – coil mat

1. Press the 'ON' button.

The illuminated display of the ON button and the polarity change display indicate that the system is ready for use.

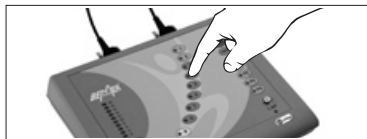2. Select the desired intensity by pressing one of the intensity levels 1 to 10 (8 minutes duration each).

- The selected intensity level will light up. Running lights for all intensity levels indicate the proper function of the system. The time display indicates the remaining therapy time in 2 minute intervals.

Instead of the intensity levels you can also select one of the four programs (P1 to P4).

P1 = 8 min., P2 = 12 min., P3 = 16 min., P4 = 20 min.

- Within the programs the intensity levels change. You can recognize the change of intensity by the indicator light of the respective intensity level. The time display shows progress in time.

At the end of the therapy period

- The running lights turn off
- An acoustic signal sounds
- After one minute without input the control unit switches off automatically

To restart the unit, repeat step 2.


If you want to stop the therapy before the programmed time, press the STOP button. If you want to disable the acoustic signal (three beeps) at the end of the treatment, keep the respective button depressed for approx. 2 seconds when activating the chosen program or level.

9.2 Control unit – intensive applicator

1. Press the 'ON' button.

The indicator light of the ON button and the polarity change display indicate that the system is ready.

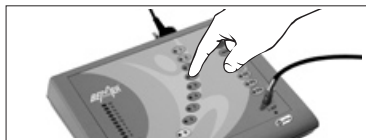2. Select the desired intensity by pressing one of the intensity levels 1 to 10 (8 minutes duration each).
 - The selected intensity level will light up. The illuminated display of the intensive applicator lights up. Running lights for all intensity levels indicate the proper function of the system. The time display indicates the remaining therapy time in 2-minute intervals.

Instead of the intensity levels you can also select one of the four programs (P1 to P4).

P1 = 8 min., P2 = 12 min., P3 = 16 min., P4 = 20 min.

- Within the programs the intensity levels change. You can recognize the change of intensity by the steadily illuminated diode of the respective intensity level. The time display shows progress in time.

At the end of the therapy period

- The running lights turn off
- An acoustic signal sounds
- The control unit switches off automatically after one minute if there is no input

To restart the unit, repeat step 2.

If you want to stop the therapy before the programmed time, press the STOP button.

If you want to disable the acoustic signal (three beeps) at the end of the treatment, keep the respective button depressed for approx. 2 seconds when activating the level or program of your choice.

An increased intensity results from the targeted treatment with the intensive applicator; however, no negative effects are known or to be expected.

9.3 BEMER Control

The BEMER Control is a handy, maintenance free, battery operated control instrument that makes the electromagnetic field created by the BEMER impulse audible. The BEMER Control can be used to check the function of the system and to demonstrate the existence and intensity of the electromagnetic field.

The contents include a 9 V block battery that is not yet connected. When the battery is discharged it can be replaced with any standard 9 V block battery.

(Please take note of Item 12.3)

1. Open the battery compartment at the back. Open the cover by gently pushing down and moving in the direction of the arrow.

When connecting the battery move the sliding switch to OFF.

2. Connect the 9 V block battery and close the battery compartment.

10. User instructions

10.1 Principles of application

The BEMER complete system features two distinct forms of application.

1. Coil mat

Therapy using the coil mat in accordance with the basic plan (BP) is aimed at a general improvement of the user's overall energy level. As a result, the conditions for a variety of metabolic processes and self-regulatory mechanisms are improved. It is the foundation for prevention and therapy.

2. Application module

- **Intensive applicator** for local application and for use when traveling (can be used with or without its pouch). Attach the pouch with the intensive applicator with the "EMF OUT" side towards the part of the body to be treated. If you do not want to use the pouch, hold the intensive applicator with the "EMF OUT" side towards the part of the body to be treated.

- **Coil cushion** (accessory) for local application and for use when traveling.

- **Multi-function applicator** (accessory) for use in therapeutic practice.

The application modules aim to achieve localized effects and are an ideal supplement to the treatment with the coil mat. Since magnetic fields decrease in intensity with increasing distance, deep effects usually require higher intensity levels.

The four intensive programs (P1 to P4) achieve a progressive 'deep treatment'.

10.2 Frequency of application

The BEMER complete system should be used daily if possible. Sporadic applications will reduce the therapeutic effect and may therefore require longer treatment. No overdose and tolerance effects are known for the BEMER complete system.

Applications can be spread out over the day as desired and depend on the purpose of the application and the time available to users.

- In general, for the basic program applications should be three to four hours apart. For preventative purposes and for maintaining good health and vitality, one to two applications per day are sufficient.
- In the physician's or therapist's practice usually only one therapy session per day is possible; as a rule, this consists of one application of the coil mat followed immediately by a specialized local application with one of the other coil applicators (double application).
- When used at home, applications can be spaced according to individual needs and schedules.
Generally, several individual applications spread throughout the day are more effective than a few double or multiple applications.
- However, sometimes schedules allow only up to two applications per day. In this case we recommend double applications where this is therapeutically required, similar to the situation in a practice.
- In case of injuries subsequent applications of the BEMER complete system may be necessary.

10.3 General instructions for the coil mat

Both sides of the coil mat and the coil cushion produce the same effect. The body should be placed in the middle of the coil mat (length 1.70 m / 5.6 feet). For taller people it does not matter whether the head, the legs, or both extend beyond the coil mat.

People with physical disabilities can be covered with the coil mat while sitting or lying down. The connection cable of the mat is attached to the top end of the coil mat (head end).

Thousands of applications in therapeutic practice have shown that people react very differently to the individual intensity levels, and frequently a correlation between metabolic activity, intensity and effect can be detected. Different reactions of circulatory parameters to different intensities of the BEMER field could be observed. The experiences gained in therapeutic practice, which are the basis for these application recommendations, are supported by photoplethysmographic studies (method for studying circulation parameters).

Experience shows that the poorer the metabolic condition of the individual, the lower the intensity level that should be used. Possible indications for a poor metabolic condition are reduced physical and / or mental capacity, fatigue, illness, pain, an increased resting pulse, shortness of breath, mental stress, and more.

10.4 Basic plan

BASIC PLAN (application using coil mat)		
	1. Application	2. Application
Week 1	Level 3	
Week 2	Level 4	
Week 3	Level 5	
Week 4	Level 6	
Week 5	Level 3	
Week 6	Level 4	
Week 7	Level 5	
Week 8	Level 6	

As of week 9, the base schedule can be repeated from weeks 5 till 8.

The basic plan aims to support regulatory mechanisms with corresponding positive physiological effects. The basic plan should be used for all conditions of illness or if well being is compromised as well as in situations when metabolic conditions are unclear.

The basic plan application (BP application) with the coil mat should take place several times a day.

It is possible to combine the coil mat with other levels or programs or with other applicators (double applications).

The basic plan includes one or several applications per day of one of the levels and one application of one of the programs. We recommend the application of one of the levels in the morning.

One cycle of BP applications lasts 4 weeks. Begin with level 3 and increase the intensity by 1 level each week, up to level 6. In addition, use program 1 once a day. In the 5th week a new cycle begins, starting again with level 3. From then on program 3 is used once a day.

Subsequent cycles correspond to the second cycle.

Note

Under certain conditions (e.g. when travelling) it is possible to use the coil cushion / the intensive applicator as part of the BP applications. For this purpose the reduction cable is used to connect the control unit with the respective coil applicator and the applicator is then positioned over the chest / stomach area (solar plexus). The coil cushion can only be used with the BEMER signal. The BEMER *plus* signal cannot be used.

In cases of extremely poor metabolic condition (e.g. strong hyperacidity), electromagnetic hypersensitivity and susceptibility to changes in weather as well as for people who are psychologically unstable, the application should start at level 1 and then increase in weekly intervals up to level 3; the basic plan can be used at a later date.

Note

Should the metabolic condition worsen during a cycle of the BP application (e.g. due to an infection, an injury / accident or an operation etc.), the current cycle should be terminated and a new cycle started using level 3 and P1.

If it is not possible to systematically increase the therapy level with the BP application (e.g. because patients can only visit their physician's or therapist's practice twice a week), the coil mat is used with P3 instead of increasing the intensity level week by week.

10.5 General effects of the different intensity levels of the coil mat

Levels 1 and 2

- Relaxes and calms the vegetative nervous system
- Reduces stress

Can be used instead of or in addition to the basic plan for the following conditions: insomnia (immediately before going to bed), severe psychological stress, depressive moods, hyperactivity, general muscle aches.

It is also suitable for: electromagnetic hypersensitivity, severe hyperacidity, nervousness and restlessness.

Levels 3 to 8

Variable range for the general support of regulatory mechanisms and/or for the improvement of various physiological metabolic processes.

For people in very good to good physical and mental condition and with good overall well being we recommend the application of the medium to high levels (6 to 8).

For most people with conditions of illness and general lack of energy we recommend levels 3 to 6 according to the basic plan.

Note: levels 3 to 4 should be used by healthy people suffering from acute short term stress situations or short term physical strain (e.g. athletic training) in the regeneration phase.

Levels 9 and 10

- Regenerating deep-tissue effect on bones and joints
- Activates and stabilizes the immune system

To be used as a supplement to the basic plan for the following conditions: osteoporosis, multiple arthroses, polyarthritis, tendency to thrombosis, auto-immune diseases, weakness of the immune system and more. It is often possible to address specific conditions with higher intensities and hence greater effectiveness using the different coil applicators (without reduction cable) on localized areas.

10.6 Additional applications

If used without the reduction cable, the coil cushion and intensive applicator achieve three times the intensity of the coil mat for the same settings of the control unit.

Since people vary a great deal in their body dimensions and, for example, the circumference of joints can vary widely due to swelling, it is not possible to recommend exact levels for the treatment of certain conditions in which a certain depth of penetration into the body would need to be achieved. For this reason, we have developed four programs in which individual intensity levels are applied in sequence and increased at certain intervals, which gradually increases the depth of penetration.

10.7 Supportive measures when using the BEMER application

- Adopt a comfortable body position, open restrictive clothing
- The BEMER applications should take place in comfortable and quiet surroundings with a comfortable temperature (approx. 22 °C / 72 degrees F, otherwise use a blanket!).
- We recommend drinking a glass of non-carbonated water before and after each BEMER application.
- As a rule, it is important to ensure sufficient hydration, i.e. approximately 2 liters (2 quarts) of non-carbonated water, spread over the day.
- No caffeinated drinks should be consumed immediately before a BEMER application.
- In addition, it is recommended that you do not smoke immediately after a BEMER application.

10.8 Programs of the BEMER complete system

PROGRAM	DURATION	DEPTH OF PENETRATION INTO ORGANISM
	8 min. 	low Treatment for superficial conditions near the skin
	12 min. 	low to medium Treatment for conditions deeper in the body
	16 min. 	
	20 min. 	deep Treatment of conditions deep in the body

Note

The indications listed in the table are only examples for certain penetration depths as the BEMER application is not a classic indication-related therapy but seeks to address the nature of illness and decreased capacity.

11. Functional test

11.1 Control unit

The control unit is functioning properly when the following light indicators are illuminated:

- The ON/OFF button
- The polarity display
- The selected intensity level
- The remaining therapy time

11.1.1 Functional test via error reporting system

The BEMER complete system is equipped with an automatic error reporting system. This is designed to function when only one applicator is connected.

Applicator failure or a bad connection are indicated visually at the control unit (running lights from levels 1 and 10 towards the center) and audibly (via a repeated signal).

11.2 Power supply unit operation

When connecting the power supply unit make sure that the operating voltage is correct: INPUT: 100 V – 240 V (do not use higher voltages). The power supply unit (type: FW7333M/15) guarantees correct function. Before connecting the power supply unit to the electrical supply it is important to ensure that the electrical supply has the same voltage as the operating voltage indicated on the power supply unit. The LED at the power supply unit will light up when properly connected to the electrical supply.

11.3 Coil mat

The control unit indicates the correct function of the coil mat.

- No error message
- Check the coil mat with the BEMER Control

11.3.1 Intensive applicator

The control unit indicates the correct function of the intensive applicator.

- No error message
- Check the intensive applicator with the BEMER Control
- Double check with optional accessories e.g. coil cushion

11.3.2 Checking function with the BEMER Control

Move the slide switch to the ON position. The control light beneath will light up. When switched on, the BEMER Control registers each electromagnetic field within 30 to 500 Hz.

Point the sensor area towards the area to be tested.

If an electromagnetic field is present, an acoustic signal will sound.

For example, the electromagnetic field around the power supply unit can be used for this functional test. If you point the BEMER Control towards the power supply unit, an acoustic signal will sound.

To turn the BEMER Control off, move the sliding switch to the OFF position. You should always turn the BEMER Control off immediately after use in order to preserve the charge of the 9 V battery.

12. Cleaning, servicing, disposal

12.1 Cleaning and disinfection

- Disconnect the power supply unit from the electrical network
- Clean the control unit, the coil mat and the intensive applicator with a slightly damp, non-abrasive cloth. Only use mild, non-abrasive cleaners in accordance with manufacturer's instructions.

None of the complete system components are waterproof or water resistant!

Do not allow moisture to penetrate.

Do not use abrasive substances or cleaners containing alkalines, acids or alcohol!

Do not use any sharp or pointed objects.

12.2 Maintenance/service

When used correctly, the BEMER complete system does not need servicing. Nevertheless, you should regularly check the condition of the equipment according to the instructions in Section 8. Start up procedures.

It is always good to know that everything functions properly; professional servicing and repair of the BEMER complete system ensures your peace of mind.

Note

For information please contact your distributor or supplier.

Servicing and repair may only be carried out by the manufacturer or an authorized agent; non-compliance will void your warranty.

Note

Before sending the BEMER complete system or some of its components to the manufacturer or an authorized agent for testing, please carry out the functional test as described in Section 10.

To ensure that the problem is not caused by an error in operation or to facilitate handling of any repair, please contact your service provider.

Use the original packaging to ensure safe shipping.

12.3 Disposal

Please dispose of the unit or batteries in accordance with local regulations. For further information please contact your supplier or the manufacturer directly.

Dispose of used batteries properly; do not throw into an open fire. Risk of explosion!
Do not dispose of batteries in regular household garbage!

13. Troubleshooting

Problems with the control unit

No indicator light

Cause	Solution
Power supply unit not connected to the electrical outlet.	Connect the power supply unit to the electrical outlet.
Power supply unit not connected to the control unit.	Connect the power supply unit to the control unit.
Control unit not switched on.	Press the ON button.
Power supply unit or control unit defective.	Send the power supply unit and / or the control unit to the manufacturer or a manufacturer authorized dealer.

No intensive applicator indicator light

Cause	Solution
Applicator is not connected.	Connect the intensive applicator to the control unit.
BEMER Control detects no signal, applicator or control unit defective.	Send the applicator and / or the control unit to the manufacturer or a manufacturer authorized dealer.

LED minute display turns off, abnormal sequence of running lights and audible signal.

Cause

No applicator connected.
(e.g. coil mat, intensive applicator)

Solution

Press the STOP button.
Connect an applicator to the control unit and select the desired treatment level.

BEMER Control detects no signal, applicator or control unit defective.

Send the applicator and /or the control unit to the manufacturer or a manufacturer authorized dealer.

Problems with the BEMER Control

No indicator light

Cause

Battery discharged.

Solution

Replace with a new 9 V block battery.

No acoustic signal

Cause

Control unit not switched on.

Solution

Press the ON button.

Connected applicator is defective.

Send the applicator to the manufacturer or a dealer authorized by the manufacturer.

BEMER Control is defective.

Send the BEMER Control to the manufacturer or a dealer authorized by the manufacturer.

14. Warranty

The manufacturer provides a warranty for the BEMER complete system covering flawless workmanship and materials under normal operating and servicing conditions for a period of 24 months from date of delivery.

This warranty does not cover recalibration and / or servicing of any kind that may be required.

The warranty only covers the first purchaser and does not cover products or components that have been used improperly or which have been subject to modifications. If the BEMER complete system is not used in accordance with its intended purpose, all warranties become void.

Liability under the warranty is limited to repairing or replacing products that have been returned to the manufacturer within the warranty period. The manufacturer must establish that the product is defective and that any such defect is not due to inappropriate handling or modifications of the unit or abnormal operating conditions.

All liabilities under the warranty become void if repairs of the BEMER complete system are carried out by persons not authorized by the manufacturer.

The above warranty provisions apply exclusively and in place of all other contractual and statutory warranty obligations, including but not limited to, the legal warranty of marketability, fitness for purpose and serviceability.

The manufacturer accepts no liability for any direct or indirect consequential damage, irrespective of whether caused by lawful, unlawful or other actions.

No liability is accepted under warranty for personal injury and damage to property where claims are based on one or several of the following causes:

- Use of the BEMER complete system not in accordance with intended purpose.
- Improper start up procedures, operation and servicing of the BEMER complete system.
- Failure to comply with notes / warnings contained in the operating instructions.
- Unauthorized modifications of the BEMER complete system.
- Failure to observe the operating instructions.
- Improper repairs and the use of spare parts or accessories not recommended by the manufacturer.
- Catastrophes and force majeure.

A claim under the warranty can only be considered if the warranty form enclosed has been completed correctly and returned to the manufacturer.

15. Technical data

Control unit BEMER Plus	Art. No. 001120
Dimensions	25.5 x 20.5 x 5 cm / 10 x 8 x 2 inches
Weight	0.75 kg / 1.65 lb
Operating voltage	12 – 15 V DC
Power input	400 mA
Power output	6 VA
Maximum level	15 V DC depending on the applicator

The output changes in accordance with the selected level
(10 levels available)

- The selected intensity is indicated by the LED
- The selected fixed program is indicated by the LED
- The remaining treatment time is indicated by the LED
- Automatic error messages
- The use of ports 1-3 is variable and can be used
for the applicator and voltage supply
- Use port I-AP front for the intensive applicator

Power supply unit (country-specific)

Power supply unit

BEMER EURO	Art. No. 001111
Primary voltage	100 – 240 V AC / 200 mA
Frequency	50 – 60 Hz

Secondary voltage	15 V DC / 530 mA
Safety rating	T40 / F
Protection rating	IP40
Protection class	SK II
Connection cable	2 m / 6.5 feet
Connection plug	
- Control unit	D-Sub
- Power supply connection	via EURO plug
Power supply indicator	via green LED

Power supply unit

BEMER int.

Art. No.	001112
Primary voltage (optional)	100 – 240 V AC / 400 mA
Frequency	50 – 60 Hz
Secondary voltage	15 V DC / 1 A
Safety rating	T40 / F
Protection rating	IP40
Protection class	SK II
Connection cable	2 m / 6.5 feet
Connection plug	
- Control unit	D-Sub
- Power supply connection	Euro plug adaptor USA plug adaptor UK plug adaptor AU plug adaptor
Power supply indicator	via green LED

Coil mat BEMER

Art. No. 001123

Dimensions

73 x 173 x 2 cm /

28.75 x 68 x 0.75 inches

Weight

2.5 kg / 5.5 lb

Surface material

100 % polyester

Coils

3 pairs of coils

Field strength (max. Ø)

35 µT

Connection cable

200 cm / 6.5 feet

Connection plug

D-Sub

Properties

elastic, suitable for rolling up

Intensive applicator BEMER Art. No. 001125

Dimensions

7 x 6 x 3 cm / 2.75 x 2.4 x 1.2 inches

Weight

200 g / 7 oz

Coil

Special coil hermetically enclosed

Field strength (max. Ø)

100 µT (EMF OUT)

Connection cable

200 cm

Connection plug

Cinch

Pouch (intensive applicator) 100 % polyester

BEMER Control

Art. No. 001302

Dimensions

6 x 9.5 x 2.5 cm /

2.4 x 3.75 x 1 inches

Weight (incl. battery)

100 g / 3.5 oz

Power supply

9 V block battery or Accu-pak

Battery life

approx. 20 hours

Power input	approx. 2 mA
Sensor range	from 1 μ T, 30–500 Hz depending on the strength of the electromagnetic field
Sound frequency	depends on the intensity of the electromagnetic field

Note

Do not store the BEMER complete system in places where it is exposed to the weather, poisonous gases or liquids.

Store the BEMER complete system in a place where it is not exposed to direct sunlight and where the specified temperature and humidity ranges are not exceeded.

16. Recommendation for optional accessories

[Art. No.: 001512]

Accu-Pak

External battery which makes mobile use of the BEMER complete system possible.

[Art. No.: 001308]

12 V car adaptor cable

Allows mobile use of the control unit and coil cushion in your vehicle.

[Art. No.: 001308]

Coil cushion

Handy applicator combining the properties of the coil mat and the intensive applicator – ideal for local application and for travel.

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